ATHLEX

**PERCUSSION
MASSAGE GUN**

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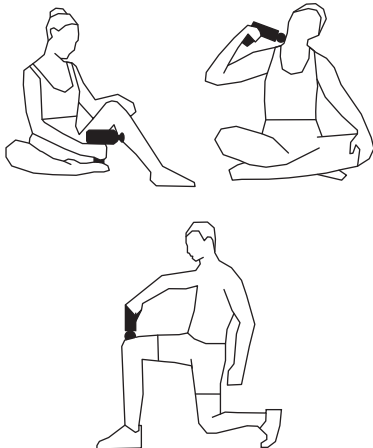
SPECIFICATIONS:

- 4 massage intensity levels
- Built-in 2,000 mAh battery
- Power: 5 VDC, 1 A
- Up to 2 hours of use per charge
- Charging time: up to 4 hours
- Amplitude: 6 mm
- Max rpm: 3,200
- Power consumption: 30 W
- 4 massage heads for various muscle groups

INCLUDED:

- Percussion massage gun
- 4 massage heads
- Charging cable
- Storage bag
- Instruction

Warranty: 2 years

The warranty does not cover wearable parts (charging cable, storage case, massage heads)

MATERIALS: ABS, aluminum

Massage heads included:


Flat

Spherical

Cylindrical

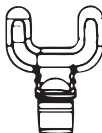U-shaped


Flat: for relaxation massage

Spherical: for large muscle groups


Cylindrical: trigger points and deep tissue massage

U-shaped: neck and Achilles tendon


How to attach massage heads:

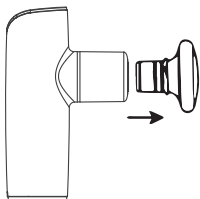remove

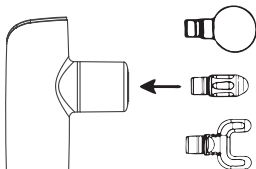install

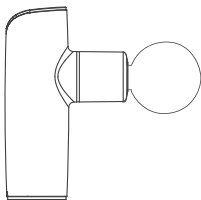installation is completed

SAFETY PRECAUTIONS:

- Read and understand the manual before using the device
- The device is intended for personal use only. Not for commercial use.
- Consult your physician before using the device
- Use this massager only as intended and as described in this manual
- Do not use the device if you are pregnant or feeling unwell
- Do not use when sleeping
- Do not use on damaged skin
- Not suitable for children or persons with disabilities
- Do not use the device in humid or dusty environments, as this may result in malfunction or electric shock
- Do not use the device at temperatures above 40 °C
- Do not expose the massager to direct sunlight or excessive heat
- Do not insert any foreign objects into the device
- Keep the device away from flammable materials
- Do not immerse the device in water.

Do not handle it with wet hands

- Do not cover the device with fabric or clothing during use, as the massager may overheat and catch fire
- Do not stand on the device, as this may result in injury or damage
- Do not disassemble the device, as this may result in injury or electric shock
- If you do not use the device for a long time, store it in a dust- and moisture-proof place
- Clean the massager with soft cloth. Do not use solvents or other aggressive cleaning agents
- In the event of a malfunction, take the device to the service center. Do not try to repair the device yourself. Failure to follow these safety precautions could void your warranty
- The manufacturer shall not be responsible for any damage caused by improper use of the device

HOW TO USE:

We recommend you fully charge the massager before the first use.

To switch on the massager, press and hold the power button for 3 seconds

To switch off the massager, also press and hold the power button for 3 seconds.

To change the mode, consecutively press the power button.

CHARGE INDICATOR:

When charging the massager, the indicator is red.


When the massager is fully charged, the indicator changes to green.

When the massager is turned on, the charge indicator flashes three times every five seconds. When the battery charge is below 60%, the indicator flashes two times every five seconds. When the battery charge is below 30%, the indicator flashes once every five seconds. When the battery charge is below 10%, the red indicator starts flashing red, and you should charge the massager.

INDICATORS:

The indicators show the current mode.

One steady light: mode 1 (1,800 rpm); two steady lights: mode 2 (2,300 rpm); three steady lights: mode 3 (2,800 rpm); three flashing lights: mode 4 (3,200 rpm).

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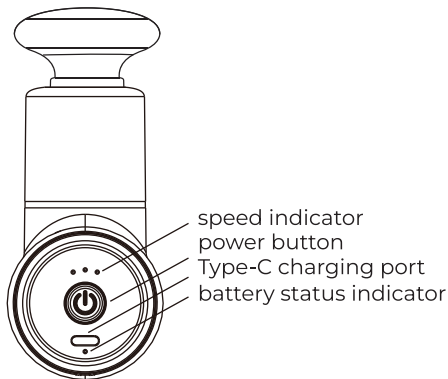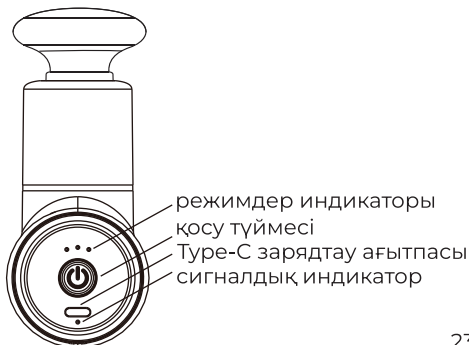
This product is legally distributed in the Russian Federation. To confirm the date of purchase for warranty repairs and for any other legal claims, please preserve the purchasing documents (check, receipt, and other documents confirming date and place of purchase). The manufacturer reserves the right to change the configuration, design, functionality of the product, its hardware and software, without prior notification.

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